

This product may be purchased from Connevans Limited secure online store at www.DeafEquipment.co.uk



DryDome.

DryDome is a compact dryer from the makers of Dry & Store®.

Gentle convection heat removes damaging moisture from hearing instruments of all types. A subtle LED indicator glows for the 8-hour cycle.

Everyday use is recommended for best results.

Features

Instant Activation: Easy-to-use 1 touch interface

High Evaporative Rate: 0.25 grams per hour

Best-in-Class 2 Year Warranty

Automatic off after 8 hours

No need to remove batteries – just open battery door

Accommodates 2 hearing aids or 1 implant processor

Physical Properties	
Weight	2.5 oz. (70.87 gm)
Unit size	3.45" (8.76 cm) diameter
Drying compartment	single tray 2.25" (5.7 cm) diameter


Electrical Characteristics

Conforms to IEC 60950-1 Ed. 2.0, CISPR 14-1, Ed. 5.1: 2009, CISPR 14-2, Ed. 1.2: 2008

Input	100-240 Volts AC, 50/60 Hz, 1 Amp
Output	5 volts DC, 1.0 Amp Max. Micro USB connector
Power supplies used with this product are certified to one or more of the following:	UL, CE, ETL, CSA, GS, PSE/JET, C-Tick, KC, and SABS. Consult label on power supply for certifications.
Operating Temperature	Up to 113°F (45°C) when using manufacturer's power supply



pry^{pc}

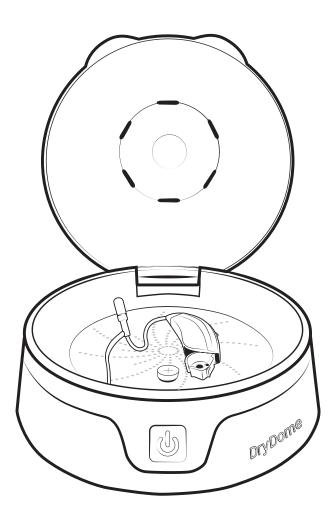

EAR TECHNOLOGY CORPORATION Helping people hear better, every day.

info@eartech.com | 1-800-327-8547 | eartech.com

P.O. BOX 1516, JOHNSON CITY, TN 37605 USA